


108900

(з вмістом 3.3 мл формаліну + тритон-Х)

Фекальний паразитологічний концентратор  
 Міні Parasер без розчинника

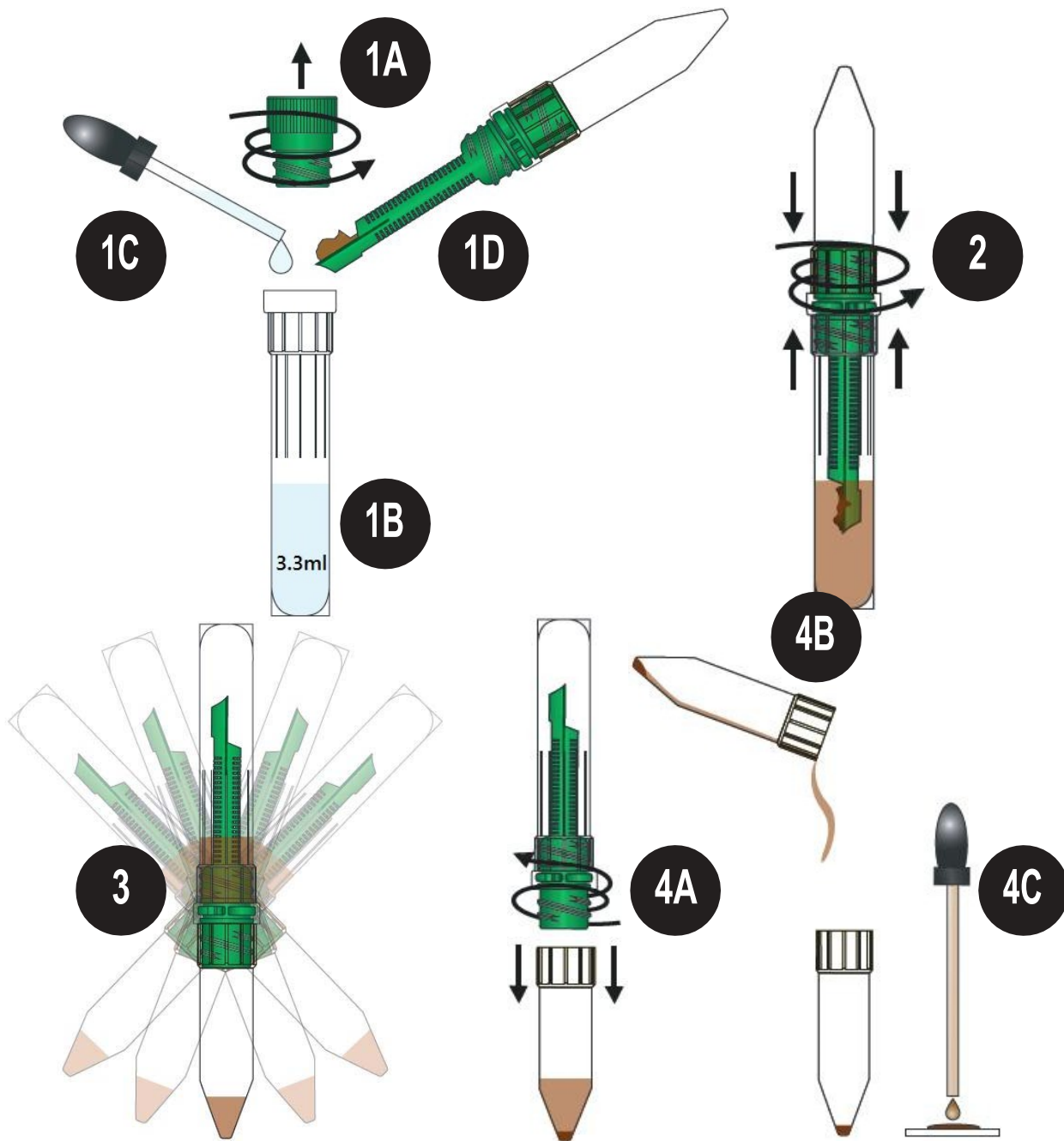

Умови зберігання та термін придатності дивитися на етикетці.

Будь ласка, дотримуйтесь наведених нижче вказівок під час роботи з Mini Parasep® SF. Щоб уникнути перехресного забруднення, концентратор Mini Parasep® SF повинен завжди бути закритим, крім випадків, коли проводиться внесення зразка або отримання концентрованого зразка для дослідження.

### **Підготовка зразка**

Якщо використовується попередньо заповнений Mini Parasep® SF, перейдіть до пункту 1D

- 1A Відкрутити кришку.
- 1B Додати 3,3 мл фіксативу.
- 1C Додати одну краплю поверхнево-активної речовини (наприклад, розчину Арасог Тритон X) у камеру.
- 1D Внести мірну ложку зі зразком калу у фіксатив за допомогою лопатки на кінці фільтру Mini Parasep® SF. Ретельно перемішати за допомогою лопатки Mini Parasep® SF. Якщо зразок твердий, розбийте його кінцем лопатки.

### **Емульгування**

- 2 Герметично закрити Mini Parasep® SF, закрутивши фільтр/конус для осадження. Вихровими рухами або струшуванням, перетворити на емульсію, тримаючи прилад конусом для осадження вгору.

### **Центрифугування**

- 3 Перевернути Mini Parasep® SF і центрифугувати при 400 g протягом 2 хвилин.  
(J. Clin. Microbiol. doi:10.1128/JCM.00838-15)

Mini Parasep® SF підходить для всіх отворів центрифуги об'ємом 15 мл.

ПРИМІТКА: Для обчислення необхідної кількості RPM (об/хв) будь-якої центрифуги

$$RPM = \sqrt{\frac{g}{1.12r}} \cdot 1000$$

RPM	Швидкість обертання ротора в об/ хв
g	Сила центрифуги (макс. 1000 g)
r	Радіус, горизонтальна відстань між верхівкою конуса для осадження та центром шпинделя, виміряна в мм

### **Тестування**

- 4A Відкрутити та викинути фільтр та трубку для змішування.
- 4B Злити всю рідину, яка утворилася над осадом.
- 4C Піпетувати одну краплю осаду на слайд і накрити покривним скельцем.  
Крім того, дотримуйтесь лабораторних стандартних правил з експлуатації для підготовки слайдів.

**10% FORMALIN SAFETY DATA SHEET**

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

**SECTION 1 IDENTIFICATION OF THE SUBSTANCE/ MIXTURE AND THE COMPANY/UNDERTAKING****1.1 Product Identifier: 10% Formalin**

145200, 145300, 145400, 145420, 145800, 145900, 1460, 146200, 146300, 146400, 108900, 108910, 148998, 149910, 151000, 900000, 903000, 905000, 908000

**1.2 Relevant identified uses of the substance or mixture and uses advised against:** laboratory chemical (in vitro diagnostic)**1.3 Details of the supplier of the Safety Data Sheet:**

Apacor Limited, Unit 5 Sapphire Centre, Fishponds Road, Wokingham, Berkshire, RG41 2QL, United Kingdom  
+44 (0) 118 979 5566  
[technical@apacor.com](mailto:technical@apacor.com)

**1.4 Emergency telephone number:**

+44 (0)118 979 5566  
(Monday-Friday 0900-1700 excluding UK Public Holidays)

**SECTION 2 HAZARDS IDENTIFICATION****2.1 Classification of the substance or mixture**

Classification according to Regulation (EC) No 1272/2008 [CLP]:

Acute toxicity, Oral (Category 4), H302

Skin sensitisation (Category 1), H317

Acute toxicity, Inhalation (Category 4), H332

Germ cell mutagenicity (Category 2), H341

Carcinogenicity (Category 1B), H350

See Section 16 for the full text of H-Statements mentioned in this Section.

**2.2 Label elements**

Labelling according to Regulation (EC) No 1272/2008 [CLP]

**Pictogram****Signal word**

**Danger**

**Hazard statement(s)**

H302 Harmful if swallowed

H317 May cause an allergic skin reaction.

H332 Harmful if inhaled.

H341 Suspected of causing genetic defects

H350 May cause cancer

Contains Formaldehyde.

**Precautionary statements:**

P280 Wear protective gloves/protective clothing/eye protection/face protection.

P301 + P312 + P330 IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell. Rinse mouth.

P308 + P313 IF exposed or concerned: Get medical advice/attention.

P333 + P313 If skin irritation or rash occurs: Get medical advice/attention.

See Section 16 for the full text of H-Statements mentioned in this Section.

**2.3 Other hazards**

This substance/mixture contains no components considered to be either persistent, bioaccumulative and toxic (PBT), or very persistent and very bioaccumulative (vPvB) at levels of 0.1% or higher.

**SECTION 3 COMPOSITION/INFORMATION ON INGREDIENTS****3.2 Mixtures****Hazardous ingredients according to Regulation (EC) No 1272/2008**

Component: **Formaldehyde**

CAS No: 50-00-0

EC No: 200-001-8

Index No: 605-001-00-5

Classification: Acute Tox. 3 (H301 + H311 + H331), Skin Corr. 1B (H314), Skin Sens. 1 (H317), Muta. 2 (H341), Carc. 1B (H350)

Concentration: < 5%

Component: **Methanol**

CAS No: 67-56-1

EC No: 200-659-6

Index No: 603-001-00-x

Registration No: 01-2119433307-44-xxxx

Classification: Flam. Liq. 2 (H225); Acute Tox. 3 (H301 + H311 + H331); STOT SE 1 H370

Concentration: < 1%

See Section 16 for the full text of H-Statements mentioned in this Section.

**SECTION 4 FIRST AID MEASURES****4.1 Description of first aid measures**

Consult a physician. Show this safety data sheet to the doctor in attendance.

**If inhaled:** If breathed in, move person into fresh air. If not breathing, give artificial respiration. Consult a physician.

**In case of skin contact:** Wash off with soap and plenty of water. Take victim immediately to hospital. Consult a physician.

**In case of eye contact:** Rinse thoroughly with plenty of water for at least 15 minutes and consult a physician.

**If swallowed:** Do NOT induce vomiting. Never give anything by mouth to an unconscious person. Rinse mouth with water. Consult a physician.

**4.2 Most important symptoms and effects, both acute and delayed**

The most important known symptoms and effects are described in the labelling (Section 2.2) and/or Section 11.

**4.3 Indication of any immediate medical attention and special treatment needed**

No data available.

**SECTION 5 FIRE FIGHTING MEASURES****5.1 Extinguishing media**

Suitable extinguishing media: Use water spray, alcohol-resistant foam, dry chemical or carbon dioxide.

## 10% FORMALIN SAFETY DATA SHEET

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

**5.2 Special hazards arising from the substance or mixture**

Carbon oxides

**5.3 Advice for firefighters**

Wear self-contained breathing apparatus and full protective gear.

**SECTION 6 ACCIDENTAL RELEASE MEASURES****6.1 Personal precautions, protective equipment and emergency procedures**

Use personal protective equipment. Avoid breathing vapours, mist or gas. Ensure adequate ventilation. Remove all sources of ignition. Evacuate personnel to safe areas. Beware of vapours accumulating to form explosive concentrations. Vapours can accumulate in low areas. For personal protection see Section 8.

**6.2 Environmental precautions**

Prevent further leakage or spillage if safe to do so. Do not let product enter drains.

**6.3 Methods and material for containment and cleaning up**

Contain spillage, and then collect and place in container for disposal according to local regulations (see Section 13). Keep in suitable, closed containers for disposal.

**6.4 Reference to other sections**

For disposal, see Section 13.

**SECTION 7 HANDLING AND STORAGE****7.1 Precautions for safe handling**

Avoid contact with skin and eyes. Avoid inhalation of vapour or mist. Keep away from sources of ignition—no smoking. Take measures to prevent the build-up of electrostatic charge. For precautions see Section 2.2.

**7.2 Conditions for safe storage, including any incompatibilities**

Store in a cool place. Keep container tightly closed in a dry and well-ventilated place. Containers which are opened must be carefully resealed and kept upright to prevent leakage.

**7.3 Specific end use(s)**

No other specific uses are specified apart from those listed in Section 1.2.

**SECTION 8 EXPOSURE CONTROLS/PERSONAL PROTECTION****8.1 Control parameters**

Exposure limits: this product, as supplied, does not contain any hazardous materials with occupational exposure limits established by the region specific regulatory bodies.

	Formaldehyde 50-00-0	Methanol 67-56-1
<b>Austria</b>	STEL: 0.5 ppm STEL: 0.6 mg/m <sup>3</sup> TWA: 0.5 ppm TWA: 0.6 mg/m <sup>3</sup>	STEL: 800 ppm STEL: 1040 mg/m <sup>3</sup> TWA: 200 ppm TWA: 260 mg/m <sup>3</sup>
<b>Belgium</b>	STEL: 0.3 ppm STEL: 0.38 mg/m <sup>3</sup>	STEL: 250 ppm STEL: 333 mg/m <sup>3</sup> TWA: 200 ppm TWA: 266 mg/m <sup>3</sup>

	Formaldehyde 50-00-0	Methanol 67-56-1
<b>Denmark</b>	STEL: 0.3 ppm STEL: 0.4 mg/m <sup>3</sup> TWA: 0.3 ppm TWA: 0.4 mg/m <sup>3</sup>	STEL: 400 ppm STEL: 520 mg/m <sup>3</sup> TWA: 200 ppm TWA: 260 mg/m <sup>3</sup>
<b>France</b>	TWA: 0.5 ppm STEL: 1 ppm	STEL: 1000 ppm STEL: 1300 mg/m <sup>3</sup> TWA: 200 ppm TWA: 260 mg/m <sup>3</sup>
<b>Germany</b>	STEL: 0.6 ppm STEL: 0.74 mg/m <sup>3</sup> TWA: 0.3 ppm TWA: 0.37 mg/m <sup>3</sup>	STEL: 800 ppm STEL: 1080 mg/m <sup>3</sup> TWA: 200 ppm TWA: 270 mg/m <sup>3</sup>
<b>Ireland</b>	STEL: 2 ppm STEL: 2.5 mg/m <sup>3</sup> TWA: 2 ppm TWA: 2.5 mg/m <sup>3</sup>	TWA: 200 ppm TWA: 260 mg/m <sup>3</sup>
<b>Italy</b>		TWA: 200 ppm TWA: 260 mg/m <sup>3</sup>
<b>Poland</b>	STEL: 1 mg/m <sup>3</sup> TWA: 0.5 mg/m <sup>3</sup>	STEL: 300 mg/m <sup>3</sup> TWA: 100 mg/m <sup>3</sup>
<b>Portugal</b>	STEL: 0.3 ppm	STEL: 250 ppm TWA: 200 ppm TWA: 260 mg/m <sup>3</sup>
<b>Spain</b>	STEL: 0.3 ppm STEL: 0.37 mg/m <sup>3</sup>	STEL: 250 ppm STEL: 333 mg/m <sup>3</sup> TWA: 200 ppm TWA: 266 mg/m <sup>3</sup>
<b>Sweden</b>	STEL: 0.6 ppm STEL: 0.74 mg/m <sup>3</sup> TWA: 0.3 ppm TWA: 0.37 mg/m <sup>3</sup>	STEL: 250 ppm STEL: 350 mg/m <sup>3</sup> TWA: 200 ppm TWA: 250 mg/m <sup>3</sup>
<b>The Netherlands</b>	STEL: 0.5 mg/m <sup>3</sup> TWA: 0.15 mg/m <sup>3</sup>	TWA: 133 mg/m <sup>3</sup>
<b>UK</b>	STEL: 2 ppm STEL: 2.5 mg/m <sup>3</sup> TWA: 2 ppm TWA: 2.5 mg/m <sup>3</sup>	STEL: 250 ppm STEL: 333 mg/m <sup>3</sup> TWA: 200 ppm TWA: 266 mg/m <sup>3</sup>

**8.2 Exposure controls****8.2.1 Appropriate engineering controls**

Handle in accordance with good industrial hygiene and safety practice. Wash hands before breaks and at the end of workday.

**8.2.2 Personal protective equipment**

(a) Eye/face protection: Tightly fitting safety goggles. Faceshield (8-inch minimum). Use equipment for eye protection tested and approved under appropriate government standards such as NIOSH (US) or EN 166(EU).

(b) Skin protection: Handle with gloves. Gloves must be inspected prior to use. Use proper glove removal technique (without touching glove's outer surface) to avoid skin contact with this product. Dispose of contaminated gloves after use in accordance with applicable laws and good laboratory practices. Wash and dry hands. The selected protective gloves should satisfy the specifications of EU Directive 89/686/EEC and the standard EN 374 derived from it.

(c) Body Protection: Complete suit protecting against chemicals. The type of protective equipment must be selected according to the concentration and amount of the dangerous substance at the specific workplace.

(d) Respiratory protection: Where risk assessment shows air-purifying respirators are appropriate use a full-face respirator with multi-purpose combination (US) or type ABEK (EN 14387) respirator cartridges as a backup to engineering controls. If the respirator is the sole means of protection, use a full-face supplied air respirator. Use respirators and components tested

**10% FORMALIN SAFETY DATA SHEET**

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

and approved under appropriate government standards such as NIOSH (US) or CEN (EU).

**8.2.3 Environmental exposure controls**

Prevent further leakage or spillage if safe to do so. Do not let product enter drains.

**SECTION 9 PHYSICAL AND CHEMICAL PROPERTIES****9.1 Information on basic physical and chemical properties**

- a) **Appearance** Form: liquid  
 b) **Odour** no data available  
 c) **Odour threshold** no data available  
 d) **pH** no data available  
 e) **Melting point / freezing point** no data available  
 f) **Initial boiling point and boiling range** 100°C at 1.013 hPa  
 g) **Flash point** 85°C  
 h) **Evaporation rate** no data available  
 i) **Flammability (solid, gas)** no data available  
 j) **Upper/lower flammability or explosive limits**  
 Upper 70% (V), Lower 7% (V)  
 k) **Vapour pressure** 53hPa at 39°C  
 l) **Vapour density** no data available  
 m) **Relative density** 1.080g/cm<sup>3</sup>  
 n) **Solubility (ies)** completely miscible  
 o) **Partition coefficient: n-octanol/water** no data available  
 p) **Auto-ignition temperature** no data available  
 q) **Decomposition temperature** no data available  
 r) **Viscosity** no data available  
 s) **Explosive properties** no data available  
 t) **Oxidising properties** no data available
- 9.2 Other information**  
 No data available.

**SECTION 10 STABILITY AND REACTIVITY****10.1 Reactivity**

No data available.

**10.2 Chemical stability**

Stable under recommended storage conditions.

**10.3 Possibility of hazardous reactions**

No data available.

**10.4 Conditions to avoid**

Heat, flames and sparks.

**10.5 Incompatible materials**

No materials to be mentioned in particular.

**10.6 Hazardous decomposition products**

Carbon oxides.

**SECTION 11 TOXICOLOGICAL INFORMATION****11.1 Information of toxicological effects**

**Acute toxicity:** no data available

**Skin corrosion/irritation:** no data available

**Serious eye damage/eye irritation:** no data available

**Respiratory or skin sensitisation:** no data available

**Germ cell mutagenicity:** no data available

**Carcinogenicity:** IARC: 1 - Group 1: Carcinogenic to humans (Formaldehyde)

**Reproductive toxicity:** no data available

**Specific target organ toxicity - single exposure:** no data available

**Specific target organ toxicity - repeated exposure:** no data available

**Aspiration hazard:** no data available

**Additional Information**

Chemical Name	
Formaldehyde	LD50 oral 600mg/kg (Rat)
	LD50 dermal 270mg/kg (Rabbit)
	LC50 inhalation 0.578mg/L (Rat) 4 h
Methanol	LD50 oral - rat - 5628mg/kg
	LC50 inhalation - rat - 4h - 83.2mg/l/4h

**SECTION 12 ECOLOGICAL INFORMATION****12.1 Toxicity**

Ecotoxicity effects: contains no substances known to be hazardous to the environment or not degradable in waste water treatment plants.

Toxicity to Fish	
Formaldehyde	0.032 - 0.226: 96 h Oncorhynchus mykiss mL/L LC50 flow-through 100- 136: 96 h Oncorhynchus mykiss mg/L LC50 static 1510: 96 h Lepomis macrochirus µg/L LC50 static 22.6 - 25.7: 96 h Pimephales promelas mg/L LC50 flow-through 23.2 - 29.7: 96 h Pimephales promelas mg/L LC50 static 41: 96 h Brachydanio rerio mg/L LC50 static
Methanol	LC50 - Pimephales promelas - 28200mg / L 96h

**Toxicity to Daphnia and other Aquatic Invertebrates**

Formaldehyde	11.3 - 18: 48 h Daphnia magna mg/L EC50 Static 2: 48 h Daphnia magna mg/L LC50
Methanol	EC50 - Daphnia magna - >10000mg/l

**12.2 Persistence and degradability**

No data available.

**12.3 Bioaccumulative potential**

No data available.

Chemical Name	log Pow
Formaldehyde	0.35
Methanol	-0.77

**12.4 Mobility in soil**

No data available.

**12.5 Results of PBT and vPvB assessment**

This substance/mixture contains no components considered to be either persistent, bioaccumulative and toxic (PBT), or very persistent and very bioaccumulative (vPvB) at levels of 0.1% or higher.

**12.6 Other adverse effects**

No data available.

**10% FORMALIN SAFETY DATA SHEET**

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

**12.7 Additional information**

None.

Amended sections are indicated by a line in the border.

The information supplied in this SDS is correct to the best of our knowledge. We do not accept any liability for loss, injury or damage, which may result from its use.

**SECTION 13 DISPOSAL CONSIDERATIONS****13.1 Waste treatment methods**

**Product:** Dispose of in accordance with all federal, state, and local regulations. This combustible material may be burned in a chemical incinerator equipped with an afterburner and scrubber. Offer surplus and non-recyclable solutions to a licensed disposal company.

**Contaminated packaging:** Dispose of as unused product.



MDSS GmbH  
Schiffaraben 41  
30175 Hanover  
Germany

**SECTION 14 TRANSPORT INFORMATION**

**IATA/DOT/IMDG/TDG:** Not regulated.

**14.1 UN number:** -

**14.2 UN proper shipping name:** -

**14.3 Transport hazard class(es):** -

**14.4 Packing group:** -

**14.5 Environmental hazards:** -

**14.6 Special precautions for user:** -

**SECTION 15 REGULATORY INFORMATION****15.1 Safety, health and environmental regulations/legislation specific for the substance or mixture**

No data available.

**15.2 Chemical Safety Assessment**

No chemical safety assessment has been carried out for this product.

**SECTION 16 OTHER INFORMATION****Full text of H-Statements referred to in Sections 2 and 3**

H225 Highly flammable liquid and vapour.

H301 + H311 + H331 Toxic if swallowed, in contact with skin or if inhaled

H302 Harmful if swallowed.

H314 Causes severe skin burns and eye damage.

H317 May cause an allergic skin reaction.

H332 Harmful if inhaled.

H341 Suspected of causing genetic defects.

H350 May cause cancer.

H370 Causes damage to organs.

Acute Tox. Acute toxicity

Carc. Carcinogenicity

Flam. Liq. Flammable liquids

Muta. Germ cell mutagenicity.

Skin Corr. Skin corrosion

Skin Sens. Skin sensitisation

STOT SE Specific target organ toxicity - single exposure

**TRITON X SOLUTION SAFETY DATA SHEET**

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

**SECTION 1 IDENTIFICATION OF THE SUBSTANCE/ MIXTURE AND THE COMPANY/UNDERTAKING****1.1 Product Identifier: 1472, 172018****TRITON X Solution**

Used at concentration of  $\leq 0.1\%$  in: 145300, 145400, 145420, 145500, 145501, 145800, 145900, 146300, 146400, 146500, 146501, 108900, 108910, 108920, 108935, 148998, 149910, 149920, 151000, 249400, 900000, 901000, 903000, 905000, 906000, 908000

**1.2 Relevant identified uses of the substance or mixture and uses advised against:** for laboratory use (in vitro diagnostic).

**1.3 Details of the supplier of the Safety Data Sheet:**

Apacor Limited, Unit 5 Sapphire Centre, Fishponds Road, Wokingham, Berkshire, RG41 2QL, United Kingdom  
+44 (0) 118 979 5566  
[technical@apacor.com](mailto:technical@apacor.com)

**1.4 Emergency telephone number:**

+44 (0)118 979 5566  
(Monday-Friday 0900-1700 excluding UK Public Holidays)

**SECTION 2 HAZARDS IDENTIFICATION****2.1 Classification of the substance or mixture**

Classification according to Regulation (EC) No 1272/2008 [CLP]:

Serious eye damage (Category1), H318

See Section 16 for the full text of H-Statements mentioned in this Section.

**2.2 Label elements**

Labelling according to Regulation (EC) No 1272/2008 [CLP]



**Pictogram**

**Signal word**

**Danger**

**Hazard statement(s)**

H318 Causes serious eye damage

**Precautionary statements:**

P280 Wear protective gloves/protective clothing/eye protection/face protection.

P305+P351+P338 IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

P313 Get medical advice/attention.

See Section 16 for the full text of H-Statements mentioned in this Section.

**2.3 Other hazards**

None known.

**SECTION 3 COMPOSITION/INFORMATION ON INGREDIENTS****3.2 Mixtures****Hazardous ingredients according to Regulation (EC) No 1272/2008**

Component: **Triton X-100** (concentration 10–20%) (included in the Candidate List of Substances of Very High Concern (SVHC) according to Regulation (EC) No 1907/2006 (REACH))

CAS No: 9002-93-1

EC No: -

A registration number is not available for this substance as the substance or its use are exempted from registration according to Article 2 REACH Regulation (EC) No 1907/2006, the annual tonnage does not require a registration or the registration is envisaged for a later registration deadline.

Classification: Acute Tox. 4 (H302); Serious Eye Dam. 1 (H318)  
Concentration: 5–10%

See Section 16 for the full text of H-Statements mentioned in this Section.

**SECTION 4 FIRST AID MEASURES****4.1 Description of first aid measures**

Consult a physician. Show this safety data sheet to the doctor in attendance.

**If inhaled:** If breathed in, move person into fresh air. If not breathing, give artificial respiration. Consult a physician.

**In case of skin contact:** Take off immediately all contaminated clothing. Rinse skin with water/shower.

**In case of eye contact:** rinse out with plenty of water. Immediately consult an ophthalmologist.

**If swallowed:** immediately make victim drink water (2 glasses at most). Consult a physician.

**4.2 Most important symptoms and effects, both acute and delayed**

Irritation and corrosion. Risk of serious damage to eyes.

**4.3 Indication of any immediate medical attention and special treatment needed**

No data available.

**SECTION 5 FIRE FIGHTING MEASURES****5.1 Extinguishing media**

Suitable extinguishing media: Use water spray, foam, dry chemical or carbon dioxide. (Use extinguishing measures that are appropriate to local circumstances and the surrounding environment.)

Unsuitable extinguishing media: For this substance/mixture no limitations of extinguishing agents are given.

**5.2 Special hazards arising from the substance or mixture**

Not combustible. Ambient fire may liberate hazardous vapours.

**5.3 Advice for firefighters**

Special protective equipment for firefighters: In the event of fire, wear self-contained breathing apparatus.

Further information: Prevent fire extinguishing water from contaminating surface water or the ground water system.

**SECTION 6 ACCIDENTAL RELEASE MEASURES****6.1 Personal precautions, protective equipment and emergency procedures**

Advice for non-emergency personnel: Do not breathe vapours, aerosols. Avoid substance contact. Ensure adequate ventilation. Evacuate the danger area, observe emergency

**TRITON X SOLUTION SAFETY DATA SHEET**

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

procedures, consult an expert.

Advice for emergency responders: Protective equipment see Section 8.

**6.2 Environmental precautions**

Prevent further leakage or spillage if safe to do so. Do not let product enter drains. Discharge into the environment must be avoided.

**6.3 Methods and material for containment and cleaning up**

Soak up with inert absorbent material and dispose of as hazardous waste. Keep in suitable, closed containers for disposal.

**6.4 Reference to other sections**

For disposal, see Section 13.

**SECTION 7 HANDLING AND STORAGE****7.1 Precautions for safe handling**

Avoid inhalation of vapour or mist. For precautions see Section 2.2.

**7.2 Conditions for safe storage, including any incompatibilities**

Tightly closed. Recommended storage temperature see product label.

**7.3 Specific end use(s)**

No other specific uses are specified apart from those listed in Section 1.2.

**SECTION 8 EXPOSURE CONTROLS/PERSONAL PROTECTION****8.1 Control parameters**

Contains no substances with occupational exposure limit values.

**8.2 Exposure controls****8.2.1 Appropriate engineering controls**

Handle in accordance with good industrial hygiene and safety practice. Wash hands before breaks and at the end of workday.

**8.2.2 Personal protective equipment**

(a) Eye/face protection: Tightly fitting safety goggles. Faceshield (8-inch minimum). Use equipment for eye protection tested and approved under appropriate government standards such as NIOSH (US) or EN 166(EU).

(b) Skin protection: Handle with gloves. Gloves must be inspected prior to use. Use proper glove removal technique (without touching glove's outer surface) to avoid skin contact with this product. Dispose of contaminated gloves after use in accordance with applicable laws and good laboratory practices. Wash and dry hands. The selected protective gloves should satisfy the specifications of EU Directive 89/686/EEC and the standard EN 374 derived from it.

(c) Body Protection: Complete suit protecting against chemicals. The type of protective equipment must be selected according to the concentration and amount of the dangerous substance at the specific workplace.

(d) Respiratory protection: Where risk assessment shows air-purifying respirators are appropriate use a full-face respirator

with multi-purpose combination (US) or type ABEK (EN 14387) respirator cartridges as a backup to engineering controls. If the respirator is the sole means of protection, use a full-face supplied air respirator. Use respirators and components tested and approved under appropriate government standards such as NIOSH (US) or CEN (EU).

**8.2.3 Environmental exposure controls**

Prevent further leakage or spillage if safe to do so. Do not let product enter drains. Discharge into the environment must be avoided.

**SECTION 9 PHYSICAL AND CHEMICAL PROPERTIES****9.1 Information on basic physical and chemical properties**

- a) **Appearance** Form: clear, liquid; Colour: light yellow
  - b) **Odour** no data available
  - c) **Odour threshold** no data available
  - d) **pH** 9.7
  - e) **Melting point / freezing point** approx. 6°C
  - f) **Initial boiling point and boiling range** 200°C
  - g) **Flash point** 251°C
  - h) **Evaporation rate** no data available
  - i) **Flammability (solid, gas)** no data available
  - j) **Upper/lower flammability or explosive limits** no data available
  - k) **Vapour pressure** <1 hPa at 25°C
  - l) **Vapour density** no data available
  - m) **Relative density** 1.070 g/cm<sup>3</sup>
  - n) **Solubility (ies)** Soluble in water
  - o) **Partition coefficient: n-octanol/water** no data available
  - p) **Auto-ignition temperature** no data available
  - q) **Decomposition temperature** no data available
  - r) **Viscosity** no data available
  - s) **Explosive properties** no data available
  - t) **Oxidising properties** no data available
- 9.2 Other information** no data available

**SECTION 10 STABILITY AND REACTIVITY****10.1 Reactivity**

No data available.

**10.2 Chemical stability**

Stable under recommended storage conditions.

**10.3 Possibility of hazardous reactions**

No data available.

**10.4 Conditions to avoid**

No data available.

**10.5 Incompatible materials**

Strong acids. Strong bases. Strong oxidizing agents.

**10.6 Hazardous decomposition products**

Other decomposition products—no data available. In the event of fire: see Section 5.



**TRITON X SOLUTION SAFETY DATA SHEET**

This Safety Datasheet complies with the requirements of Regulation (EC) No 1907/2006

**SECTION 11 TOXICOLOGICAL INFORMATION****11.1 Information of toxicological effects****Acute toxicity:** no data available**Skin corrosion/irritation:** no data available**Serious eye damage/eye irritation:** no data available**Respiratory or skin sensitisation:** no data available**Germ cell mutagenicity:** no data available**Carcinogenicity:** IARC: no component of this product present at levels greater than 0.1% is identified as probable, possible or confirmed human carcinogen by IARC.**Reproductive toxicity:** no data available**Specific target organ toxicity - single exposure:** no data available**Specific target organ toxicity - repeated exposure:** no data available**Aspiration hazard:** no data available**Additional information:** RTECS: not available. To the best of our knowledge, the chemical, physical, and toxicological properties have not been thoroughly investigated.**11.2 Further information**

Triton X-100

Acute oral toxicity: LD50 Rat: 1,800 mg/kg (RTECS)

Germ cell mutagenicity: Genotoxicity in vitro Mutagenicity (mammal cell test): Mouse lymphoma test Result: negative

**SECTION 12 ECOLOGICAL INFORMATION****12.1 Toxicity**

No data available.

**12.2 Persistence and degradability**

No data available.

**12.3 Bioaccumulative potential**

No data available.

**12.4 Mobility in soil**

No data available.

**12.5 Results of PBT and vPvB assessment**

PBT/vPvB assessment not available as chemical safety assessment not required/not conducted.

**12.6 Other adverse effects**

Discharge into the environment must be avoided.

Components: Triton X-100

Toxicity to fish

LC50 Lepomis macrochirus: 2,800 - 3,200 µg/l; 96 h

Toxicity to daphnia and other aquatic invertebrates

LC50 Daphnia magna: 11.2 mg/l; 48 h

**12.7 Additional information**

No data available.

**SECTION 13 DISPOSAL CONSIDERATIONS****13.1 Waste treatment methods****Product:** Offer surplus and non-recyclable solutions to a licensed disposal company.**Contaminated packaging:** Dispose of as unused product.**SECTION 14 TRANSPORT INFORMATION****IATA/DOT/IMDG/TDG:** Not regulated.**14.1 UN number:** -**14.2 UN proper shipping name:** -**14.3 Transport hazard class(es):** -**14.4 Packing group:** -**14.5 Environmental hazards:** -**14.6 Special precautions for user:** -**SECTION 15 REGULATORY INFORMATION****15.1 Safety, health and environmental regulations/legislation specific for the substance or mixture**

Substances of very high concern (SVHC): This product does contain substances of very high concern above the respective regulatory limit (&gt;0.1% w/w), Regulation (EC) No 1907/2006 (REACH), Article 57).

Contains: Triton X-100.

**15.2 Chemical Safety Assessment**

No chemical safety assessment has been carried out for this product.

**SECTION 16 OTHER INFORMATION****Full text of H-Statements referred to in Sections 2 and 3**

H302 Harmful if swallowed

H318 Causes serious eye damage

Acute Tox. Acute Toxicity

Serious Eye Dam. Serious Eye Damage

Amended sections are indicated by a line in the border.

The information supplied in this SDS is correct to the best of our knowledge. We do not accept any liability for loss, injury or damage, which may result from its use.



MDSS GmbH  
Schiffaraben 41  
30175 Hanover  
Germany

Фекальний паразитологічний концентратор  
MiniParasep® без розчинника



Кат. номер	Продукт
108900	Mini Parasep® SF Formalin & Triton X

**Утилізувати клінічні відходи відповідно до стандартних та місцевих процедур експлуатації.**

Продукцію можна замовити безпосередньо у Аракор або у призначеного дистриб'ютора  
Щоб отримати найсвіжішу інформацію, відвідайте наш веб-сайт [www.apacor.com](http://www.apacor.com) або  
електронну пошту: [orders@apacor.com](mailto:orders@apacor.com)



ВИРОБНИК  
АПАКОР ЛІМІТЕД  
ЮНІТ 5 САПФІР ЦЕНТР  
ФІШПОУНДС РОУД, ВОКІНГЕМ  
БЕРКШИР, RG41 2QL, АНГЛІЯ  
ТЕЛ.: +44 (0)118 979 5566  
ФАКС: +44 (0)118 979 5186



**УПОВНОВАЖЕНИЙ ПРЕДСТАВНИК**

ТОВ «ДІАМЕБ ТРЕЙД»  
вул. Симона Петлюри, 25  
м. Івано-Франківськ, 76014  
тел.: +38 (0342) 775 122  
факс: +38 (0342) 775 123  
e-пошта: [info@diameb.ua](mailto:info@diameb.ua)  
[www.diameb.ua](http://www.diameb.ua)